


Laminar Flow Cabinet

Laminar Flow Cabinet

Laminar Air flow Clean bench can supply you the high cleanliness and germfree laminar air flow surroundings. All the operation are under constant positive pressure in order to prevent the infiltration of contaminated room air. The air flow can be horizontal flow and vertical laminar flow.

Air flow type

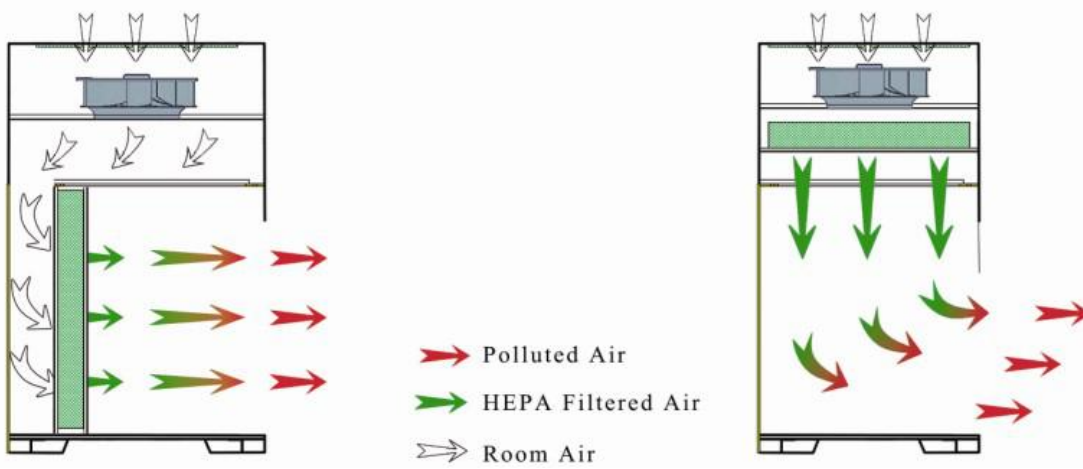

Horizontal Type

Vertical Type



BBS-H1300/BBS-H1800



BBS-H1100/BBS-H1500

Main Technical Parameters for Horizontal Type

Model	BBS-H1100	BBS-H1500	BBS-H1300	BBS-H1800
External Size (W x D x H)	1060*800*1690	1460*800*1690	1300*825*2000	1800*825*2010
Internal Size (W x D x H)	980*500*600	1380*500*600	1200*500*570	1700*500*570
Work Surface Height	700mm		720	
Display	LED		LCD	
Material	Main Body	Cold-rolled steel with anti-bacteria powder coating		
	Work Table	304 stainless steel		
	Side And Front Window	5 mm toughened glass, anti-ultraviolet radiation		
Pre-Filter	Polyester fiber, Washable			
HEPA Filter	Efficiency ≥ 99.999 % at 0.3 μ m			
Airflow Velocity	Average of 0.3~0.5m/s			
Noise	< 60 dB			
Fluorescent Lamp	21W*1	28W*1	28W*1	28W*1
UV Lamp	20W*1	30W*1	30W*1	30W*1
Consumption	320W	620W	400W	700W
Caster	Universal Wheel with leveling feet			
Power Supply	110/220V ± 10%; 50/60HZ			
Standard Accessory	Fluorescent Lamp,UV Lamp*2,Base Stand,Waterproof Socket*2			